

# Min-Max™

SKU #7046

**Low-cost temperature monitoring.**

- Low-cost temperature thermometers
- Refrigerator/freezer thermometer.
- Digital temperature range -35° C to +85° C / -35° F to +185° F
- Factory set to either Fahrenheit or Celsius.
- Press start button to refresh display.
- Min/Max reset. Freeze display and reset.
- Stand or hang on wall.
- 365 day battery life.

## Temperature Range

-35°C to +85°C /  
-31°F to +185°F

## Temperature Accuracy

(NIST-traceable)  
±0.3°C (±.54°F)  
from 0° to 10° C  
(32° to 50° F);  
±1.0°C (±1.8°F)  
at extremes

## Size:

Length: 5.5 cm/2.2 in.  
Width: 2.7cm/1.06 in.  
Depth: 1.4cm/0.55 in.

**Weight:** 32 g / 1.13 oz.

# Min-Max Plus™

SKU #7047

**Compact, economical temperature monitoring with built-in external probe and real-time LCD display.**



**The Min-Max Plus** is a reliable cost-effective digital temperature recorder that is ideal for laboratories, blood banks needing low-cost real-time information. It has an external 2 meter (6.56 ft.) probe that can measure temperatures as low as -35°C

**LCD Display:** View ambient or probed temperatures in real-time, prior to downloading data.

**Freeze Display:** To freeze the display simply press the START button. To reset: simply press once again to reset and resume.

**Visual Alarms:** This data logger features user programmable alarm limits. Should the temperature reach levels outside these limits, the LCD will display an alarm. Alarms are visible until reset by the user. Operation is as easy as pressing the START button on the logger.

**USB Interface:** Put the Min-Max Plus into the reader station to download recorded information.

**Small Form Factor:** The thumb-sized format fits easily into small spaces.


**Designed with Continuous Looping Memory:** rolls over every 30 days

**Start Only:** No stop function.

**MARATHON**  
 **PRODUCTS, INC.**

**Don't ship without us!®**

**MaxiThermal2 Software Highlights:** Easy-to-use software provides clear graphs. Report Summary & Statistics File Manager, Note Attachment, Zoom, Elapsed Time or Date/Time View, Recorded View, Celsius or Fahrenheit. Data can be exported to other Windows programs such as spreadsheets, word-processing or databases. Available for Windows 7, 8 and 10. An Enterprise Software version and 21CFR Part 11 Compliant Software are also available.

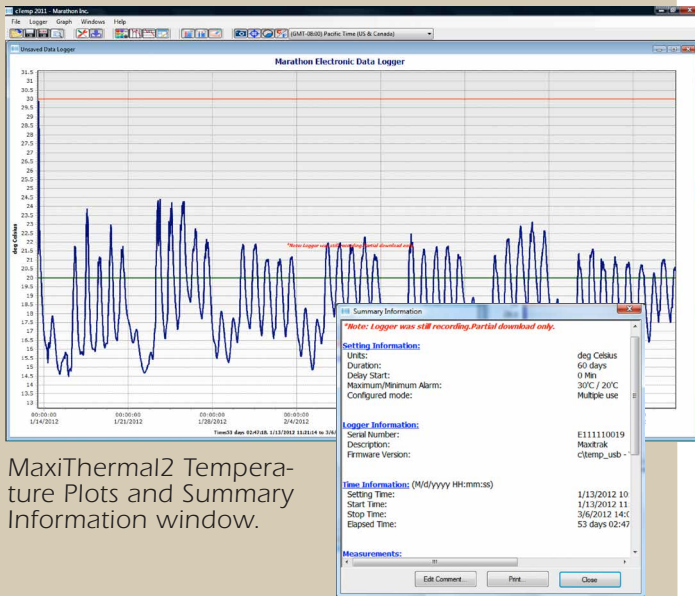

## Min-Max Plus Specifications SKU#7047

- Memory EEPROM, Storage: 8K (approx. 8,000 measurements)
- Temperature Range: -35°C to 85°C / -31°F to 185°F
- Temperature Accuracy:  $\pm 0.3^{\circ}\text{C}$  at  $0^{\circ}\text{C}$  to  $10^{\circ}\text{C}$  /  $\pm 1^{\circ}\text{C}$  at extremes  
 $\pm 0.54^{\circ}\text{F}$  at  $32^{\circ}\text{F}$  to  $50^{\circ}\text{F}$  /  $\pm 1.8^{\circ}\text{F}$  at extremes
- Measurement intervals: User programmable from 2 seconds to months
- Multi-use Temperature sensor: NTC probe with wired glass thermistor, SST tube dia 3mm 2cm long-IP67
- Probe: Silicone cable dia 2.0mm(0.08 in), 2.0 meter ( 6.56 feet) long
- Indicators: Start/Stop and Alarm
- High and low alarm limits: LCD display, ALARM set by temperature and time
- Computer interface: USB Standard Plug A to Mini B cable and Reader Station
- Internal firmware provides: Times, Start time, Interval, User description, Alarm values, Recording status, Measurement data, Firmware version, Serial number, Model type
- Power source: Lithium battery ( 1 year operating, 3 years shelf life )
- Size: Length                      Width                      Depth  
5.5 cm/2.2 in.    2.7 cm/1.06 in.    1.4 cm/0.55 in.
- Weight: 32 g / 1.13 oz with probe
- Required Accessories: MaxiThermal2 software  
Maxilog MaxiThermal2 Reader Station



**Maxilog MaxiThermal2  
Reader Station**

**Please call 1-510-562-6450**  
for more information or to place an order.

**www.marathonproducts.com**

14500 Doolittle Drive, San Leandro, CA 94577 ▪ jnakagawa@marathonproducts.com

Specifications subject to change without prior notice.